


# Staytex<sup>TM</sup>

Elastic Tubular Dressing

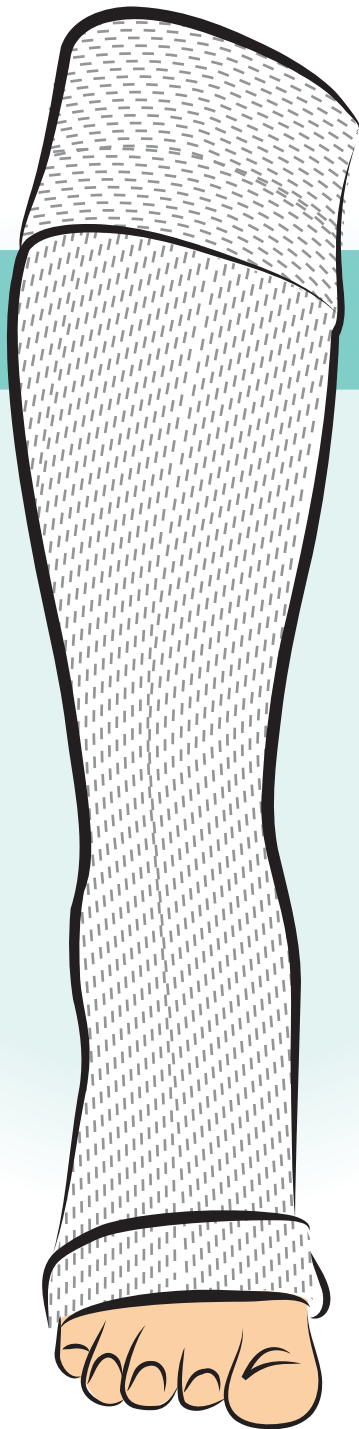

## Features:

- Secures wound dressing to affected site with light compression.
- Maintains the wound dressing in place until next wound dressing change.
- Easy to apply and reapply
- Comfortable
- Washable and Reuseable
- Eliminates the need for cumbersome and ineffective wraps.
- Anasept<sup>®</sup> Antimicrobial Skin & Wound Cleanser and Anasept<sup>®</sup> Antimicrobial Skin & Wound Gel can be used in combination with Staytex<sup>TM</sup> Elastic Tubular Dressing to maximized the properties of both products.
- Latex Free




## Elastic Tubular Dressing

### DIRECTIONS FOR USE:


- Cleanse the open wound with Anasept® Antimicrobial Skin and Wound Cleanser



- Apply Anasept® Antimicrobial Skin and Wound Gel to wound bed.

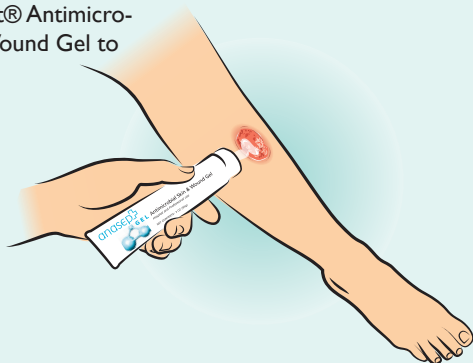

- Cover wound with sterile dressing

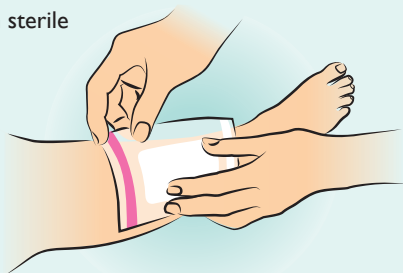


- Roll the Staytex™ tubular dressing upon your hands like you would a stocking before putting it on.

- Cut to needed length




- Place Staytex™ over the foot until toes are exposed then unroll the rolled end up the length of the leg.

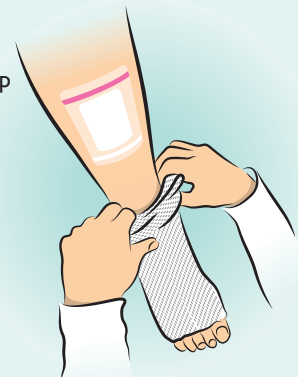

- To prevent Staytex™ tubular dressing from rolling down or up on the limb, turn down the top of Staytex™ making a three inch cuff.

- You can also turn Staytex™ inside out to prevent rolling.

- Remove daily to change wound dressing and reposition Staytex™

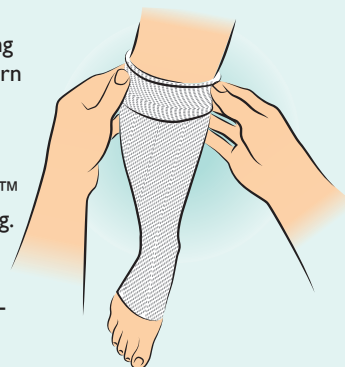

Anasept® and Staytex™ are trademarks of Anacapa Technologies, Inc.®  
Anasept® is manufactured by  
Anacapa Technologies®  
San Dimas, California  
Orders/Information: **800-489-2591**  
[www.anacapa-tech.net](http://www.anacapa-tech.net)

Staytex™ is manufactured by  
Compression Dynamics, LLC®  
Omaha, Nebraska  
Orders/Information: (402) 660-9834  
[www.CompressionDynamics.com](http://www.CompressionDynamics.com)

### RECOMMENDED CARE:

- Machine wash Staytex™ in cold water.
- Hang dry
- Remove bloodstains with Hydrogen Peroxide.